



**THOMAS
MEDICAL**
a MedGyn company



Product Catalog

Leading the Way for
Quality Infertility Products

Helping Families Grow Since 1984

Distributed in New Zealand by:



Cooper Medical Limited
57 Mahunga Drive
Mangere Bridge
Auckland 2022

www.coopermedical.co.nz - admin@coopermedical.co.nz - Ph: 0800 800 823



Who We Are

Thomas Medical, a MedGyn Company, manufactures single-use medical devices and accessories for women. We are dedicated to producing infertility diagnostic products for healthcare professionals. Thomas Medical's high-quality family of products include HS/HSG Catheters and Procedure Trays, IUI Catheters, Uterine Manipulators/Injectors, LEEP/LLETZ Electrodes, and more.



Our Mission and Vision

As a dedicated partner for medical professionals, our organization is committed to providing exceptional infertility and diagnostic products for women's healthcare. With compassion and competence, we collaborate with providers to help families grow.



What Sets Us Apart

Made in the U.S.A, Thomas Medical products meet the highest quality standards and are individually hand-inspected so healthcare practitioners can be confident that they have a reliable product for Hysterosalpingography (HSG), Saline Infusion Sonography (SIS), Sonohysterography (SHG), and additional essential diagnostic and laparoscopic procedures.

Proudly Made in the U.S.A.



Made in the USA

ISO 13485 and FDA Registered

Our Products Do Not Contain Natural Rubber Latex, DEHP or PVC

DIAGNOSTIC GYNECOLOGY

Flexible & Shapeable HS Balloon Catheters	3
Miller Advance HS Balloon Catheter	4
Steen Open-Tip Catheter	4
HSG Procedure Tray	6 - 7
eZinject Sonohysterography Catheter	8

SURGICAL GYNECOLOGY

Kronner Uterine Manipulator/Injector	9
Uterine Manipulator/Injector	10
LEEP/LLETZ Ball Electrodes	11-12
LEEP/LLETZ Square Electrodes	11-12
LEEP/LLETZ Loop Electrodes	11-12


ADDITIONAL CARE

Intrauterine Insemination Catheter	13
Disposable Uterine Sound	14
Disposable Os Finder	15
Endotracheal Tube Holder	16



HS BALLOON CATHETER

DIAGNOSTIC GYNECOLOGY

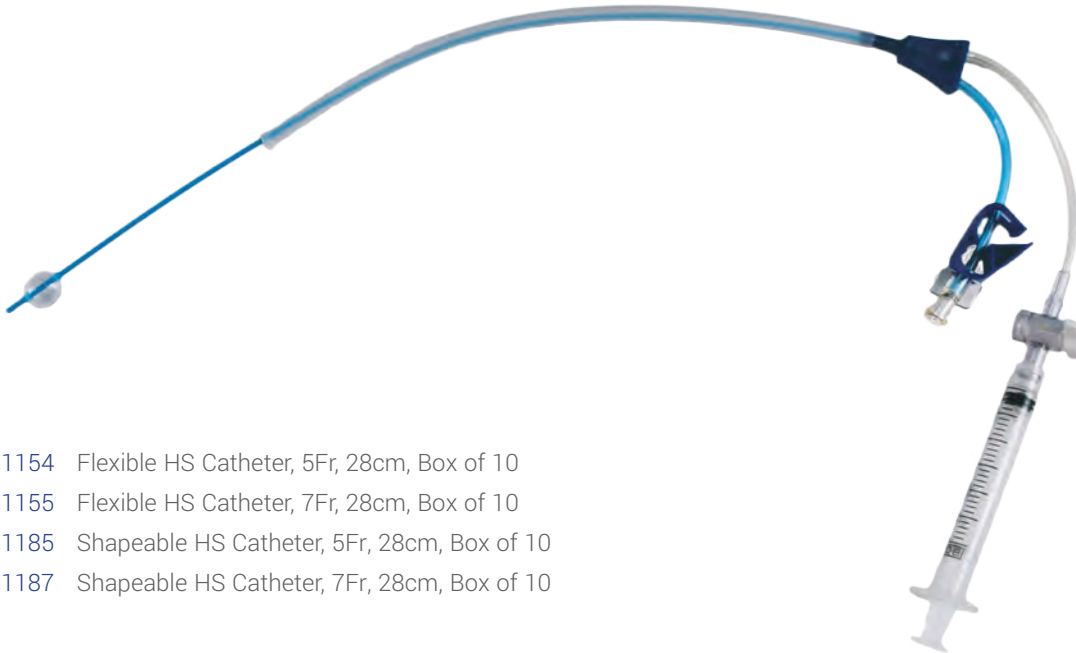

- TMI 1154 Flexible HS Catheter, 5Fr, 28cm, Box of 10
- TMI 1155 Flexible HS Catheter, 7Fr, 28cm, Box of 10
- TMI 1185 Shapeable HS Catheter, 5Fr, 28cm, Box of 10
- TMI 1187 Shapeable HS Catheter, 7Fr, 28cm, Box of 10

FLEXIBLE & SHAPEABLE HS CATHETERS

Thomas Medical's HS Catheters are designed with a smooth tip to reduce the risk of trauma upon insertion and removal. The malleable placement sheath provides stability and decreases the need for a tenaculum. A soft, concentric balloon seals the internal os and minimizes leakage of sterile saline or water-based contrast media. Further, the Shapeable HS Catheter is equipped with an embedded stylet for cases of stenosis or where the vaginal tract and cervical canal are not aligned. The shapeable shaft also provides additional rigidity.

HS BALLOON CATHETER


DIAGNOSTIC GYNECOLOGY



MILLER ADVANCE HS CATHETER

The Miller Advance HS Catheter is designed to improve patient comfort with a soft catheter tip and slender shaft that reduces the need for dilation. The sheath is shapeable to ensure maximum comfort and effectiveness.

TMI 2002 Miller Advance Catheter, 5Fr, 28cm, Box of 10

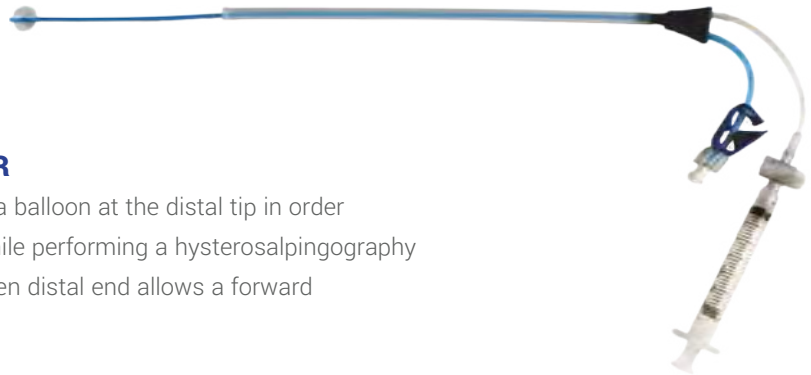

STEEN OPEN-TIP CATHETER

The Steen Open-Tip HS Catheter utilizes a balloon at the distal tip in order to provide a clearer view of the uterus while performing a hysterosalpingography or sonohysterography procedure. The open distal end allows a forward injection into the uterus.

TMI 1188 Steen Open-Tip Catheter, 5Fr, 28cm, Box of 10

HS BALLOON CATHETER

DIAGNOSTIC GYNECOLOGY

Catheter Description	Catheter Shaft	Sheath Type	Part Number	French Size	Catheter Length	Syringe Size / Balloon Volume
HS Catheter	Flexible	Shapeable	TMI1154	5Fr	28cm	3ml/1.5ml
	Flexible	Shapeable	TMI1155	7Fr	28cm	5ml/3ml
Shapeable HS Catheter	Shapeable	Shapeable	TMI1185	5Fr	28cm	3ml/1.5ml
	Shapeable	Shapeable	TMI1187	7Fr	28cm	5ml/3ml
Miller Advance Catheter	Flexible	Shapeable	TMI2002	5Fr	28cm	3ml/1.5ml
Steen Open-Tip Catheter	Flexible	Shapeable	TMI1188	5Fr	28cm	3ml/1.5ml

Each HS Catheter is Individually Hand-Inspected

Sterilization Method for HS Catheters: GAMMA



**DISTAL
TIP DESIGN**

MILLER ADVANCE
MALLEABLE
THIN SHAFT

STEEN OPEN-TIP
MINIMAL DISTAL TIP

HS CATHETER

HSG PROCEDURE TRAY

DIAGNOSTIC GYNECOLOGY

TRAY CONTENTS

1. Extension tubing, .050" ID x 36" M/F LL	1 ea.	10. Pad, Maxi	1 ea.
2. Syringe, 20CC L/L	1 ea.	11. Swab, Cotton 8"	3 ea.
3. PVP Solution 1oz	2 ea.	12. Bag, Header 10" x 24"	1 ea.
4. Jelly, Lubricating 3g.	1 pk.	13. MedGyn Disposable LED Speculum	1 ea.
5. Ointment, PVP 1g	1 ea.	14. Wrap, Poly 30" x 30"	1 ea.
6. Drape with 2 1/2 Fenestration, 18" x 26"	1 ea.	15. Tray, 16" x 6.6" x 2.3"	1 ea.
7. Needle, 18 G x 1 1/2	1 ea.	16. Tray, 5.3" x 4.5" x 2"	1 ea.
8. Gauze, 3"x 3" 4ply	5 ea.	17. Cervical Dilator	1 ea.
9. Sterile Gauze, 3x3 12ply	2 ea.	18. Sponge Stick, 6"	2 ea.

PREP IN ONE STERILE FIELD

Thomas Medical HSG Procedure Trays contain PVP solution, lubricating jelly, and PVP ointment, reducing procedure prep and set up time.

MEDGYN DISPOSABLE LED SPECULUM

MedGyn's LED Speculums are made from clear polystyrene for easy, one-hand operation and insertion. The bright, white light effortlessly improves visualization of the exam area.

PVC FREE CATHETERS

All Thomas Medical HS Catheters are PVC free and can be purchased individually or included as part of the HSG Procedure Tray.

Part Number	Procedure Tray Description	Catheter Description	Box Qty	Sterilization Method
TM1157	HSG Procedure Tray without Catheter	No catheter	10 ea.	EO
TM1160	HSG Procedure Tray with TMI1154 Catheter	5Fr HS Balloon Catheter	10 ea.	EO
TM1161	HSG Procedure Tray with TMI1155 Catheter	7Fr HS Balloon Catheter	10 ea.	EO
TM1198	HSG Procedure Tray with TMI1185 Catheter	5Fr Shapeable HS Balloon Catheter	10 ea.	EO
TM1199	HSG Procedure Tray with TMI1187 Catheter	7Fr Shapeable HS Balloon Catheter	10 ea.	EO

EZINJECT SONOHYSTEROGRAPHY CATHETER

DIAGNOSTIC GYNECOLOGY



TMI 1162-HF eZinject Catheter, Box of 10


FEATURES

The Thomas Medical eZinject Catheter is compatible with oil-based contrast media. The catheter's tip is specifically designed to prevent leakage. Rigid cervical flange and sheath for easy placement and single operator use. Long tip at distal end improves placement in cervical canal. The wide tip opening is designed for increased flow. Balloonless catheter allows for a more complete view of the uterine cavity during procedures.

Catheter Description	Shaft Type	Part Number	Tip O.D.	Catheter Length	Box Qty	Sterilization Method
eZinject Catheter	Rigid	TMI1162-HF	2mm to 3.7mm	30cm	10 ea.	EO

UTERINE MANIPULATOR/INJECTOR

SURGICAL GYNECOLOGY

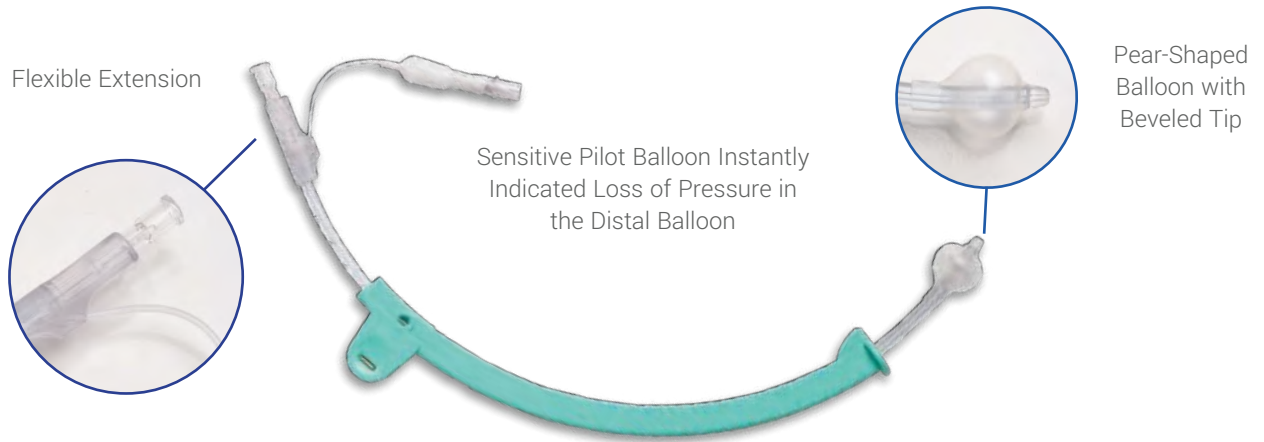

FEATURES

This unique, pear-shaped balloon conforms to the uterine anatomy and secures the internal cervical os to improve manipulation. Flexible extension reduces the risk of cracking the shaft by providing strain relief. The atraumatic beveled catheter tip reduces trauma to the uterine wall during manipulation and the beveled tip reduces trauma to the cervical canal during insertion.

KRONNER MANIPULATOR/INJECTOR (K-UMI)

The pinch-clamp on the handle allows for quick and effective repositioning during procedures.

TMI 6003 K-UMI Uterine Manipulator/Injector, Box of 12



UTERINE MANIPULATOR/INJECTOR

SURGICAL GYNECOLOGY



UTERINE MANIPULATOR/INJECTOR

The snap-on handle allows the practitioner to easily place the cervical stop against the cervical face.




TMI 1151 Uterine Manipulator/Injector, Box of 12

Description	Shaft Type	Part Number	Shaft Diameter	Handle Style	Box Qty	Sterilization Method
K-UMI Manipulator/Injector	Rigid	TMI6003	5.0mm	Pinch Clamp	12 ea.	EO
Uterine Manipulator/Injector	Rigid	TMI1151	4.5mm	Snap-On	12 ea.	EO

LEEP/LLETZ ELECTRODE

SURGICAL GYNECOLOGY

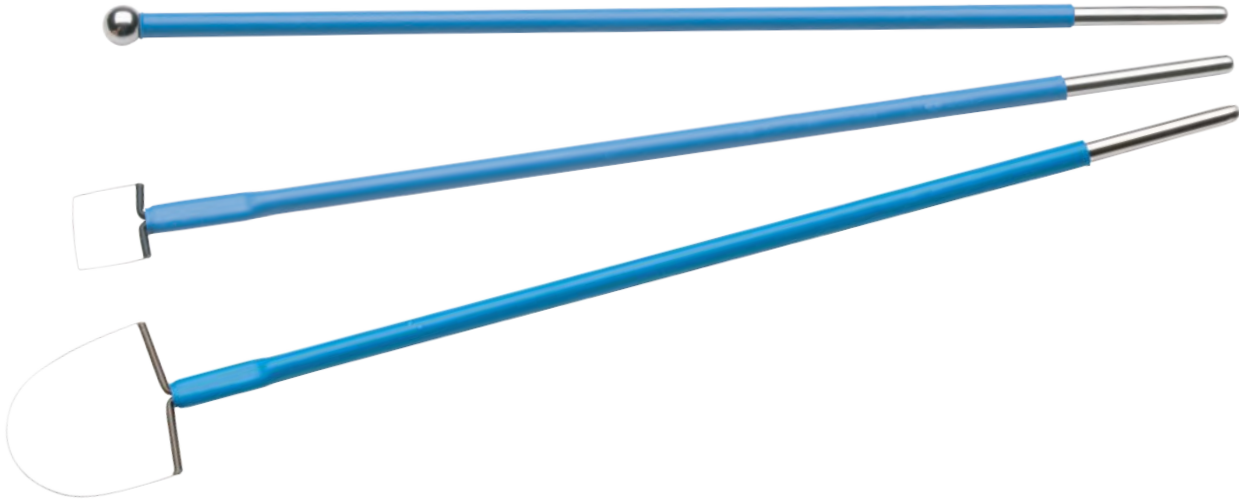

BALL ELECTRODES

- LE-13-003 3mm Diameter Ball Electrode
- LE-14-005 5mm Diameter Ball Electrode

SQUARE ELECTRODES

- LE-06-108 10mm x 8mm, Square Electrode
- LE-06-101 10mm x 10mm, Square Electrode

LOOP ELECTRODES

- LE-04-105 10mm x 5mm Loop Electrode
- LE-05-101 10mm x 10mm Loop Electrode
- LE-07-158 15mm x 8mm Loop Electrode
- LE-26-150 15mm x 10mm Loop Electrode
- LE-09-208 20mm x 8mm Loop Electrode
- LE-24-201 20mm x 10mm Loop Electrode
- LE-11-205 20mm x 15mm Loop Electrode



LEEP/LETZ ELECTRODE

SURGICAL GYNECOLOGY

FEATURES

11 different sizes ranging from
5 x 5mm to 25 x 10mm Loop, Ball and
Square Electrodes

Shaft Length 13cm (5")
Sterilization Method EO
Box Quantity 20 ea.

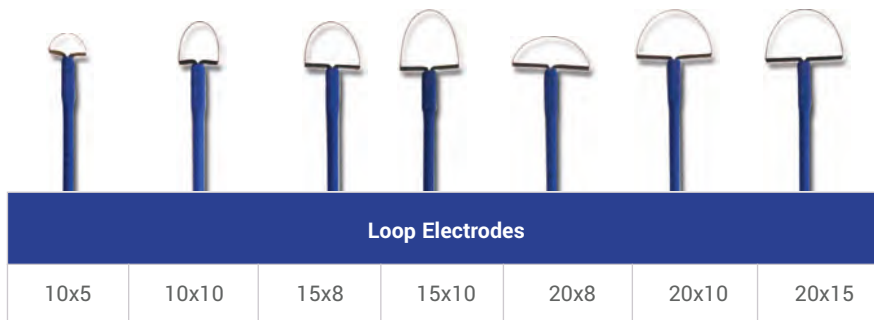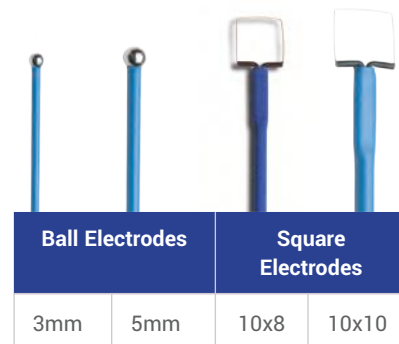

INTRAUTERINE INSEMINATION CATHETER

ADDITIONAL CARE



FEATURES

The Intrauterine Insemination Catheter is a sterile, disposable, flexible catheter designed for use in intrauterine insemination. The soft material and smooth, rounded tip help avoid damage to the uterine wall, yet the rigidity of the stylet makes entry through the cervical canal simple. The diameter of the Mini and Mini Curved helps minimize uterine cramping and maximizes sperm use while leaving minimal residue.

MINI




Mini

O.D. 1.55mm +/- .2

I.D. .8mm +/- .2

MINI CURVED

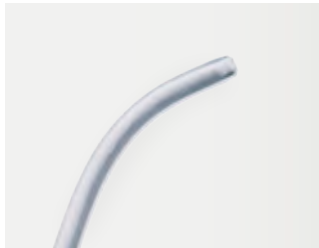

Mini Curved

O.D. 1.55mm +/- .2

I.D. .8mm +/- .2

REGULAR



Regular

O.D. 2.13mm +/- .2

I.D. 1.43mm +/- .2

022724 Mini Intrauterine Insemination Catheter, Box of 25

022725 Mini Curved Intrauterine Insemination Catheter, Box of 25

022723 Regular Straight Intrauterine Insemination Catheter, Box of 25

DISPOSABLE UTERINE SOUND

ADDITIONAL CARE

DISPOSABLE UTERINE SOUND

Made in the USA, Thomas Medical's uterine sound is simple and straight to the point. Each uterine sound is individually checked with the highest quality standard in the industry, delivering a reliable experience.

FEATURES

Sterile, Disposable, Ready to Use

Length: 24.8 cm (9 3/4")

OD Size: 3.5mm (1/10")

Shapeable Tip

DEHP Free

This Product Does Not Contain Latex

US-01-025 Disposable Uterine Sound, Box of 25



DISPOSABLE OS FINDER

ADDITIONAL CARE

DISPOSABLE OS FINDER

Gently opens the cervical os. Ideal for cervical dilation, uterine sounding or other intrauterine procedures.

FEATURES

Streile, Disposable, Ready to Use

Length: 21.5cm

OD Size: 2mm - 4mm Taper

DEHP Free

This Product Does Not Contain Latex



OS-01-050 Disposable Os Finder, Box of 50

ENDOTRACHEAL TUBE HOLDER

ADDITIONAL CARE



ENDOTRACHEAL TUBE HOLDER

Disposable and ready to use

Does not contain natural rubber latex

No tape is required to place the ET Tube Holder

May eliminate need for a bite block

Wide cotton-lined neckband minimizes skin irritation, providing maximum patient comfort

FEATURES

Endotracheal Tube Holder is used with ventilators and other breathing assistance devices.

TMI-ET1 Endotracheal Tube Holder, Box of 25

CORPORATE AFFILIATE

THOMAS MEDICAL JOINED MEDGYN

In the fall of 2020, Thomas Medical, Inc. joined with MedGyn Products, Inc., a world leader in women's healthcare products.

Because of this partnership, Thomas Medical can serve our customers and their patients better than ever. The range of clinics and communities we can reach will increase exponentially, as will our portfolio of products and services.

While Thomas Medical continues to lead the way in solutions to gynecologic diagnostic imaging, surgical gynecology, and infertility, MedGyn brings access to therapies such as family planning and gynecological oncology.

Thomas Medical is proud to join forces with this respected company to serve the women's healthcare community.





Proudly Made in the U.S.A.

100 W. Industrial Road
Addison, IL 60101 USA
T: (800) 556.0349
F: (800) 411.3754
sales@thomasmedical.com
thomasmedical.com

Follow us on Social Media



THMEZCAT-PF1



Made in the USA
2797 ISO 13485 and FDA Registered