



MASTEL SURGICAL

Where Quality Meets Precision

Distributed in New Zealand by:



Cooper Medical Limited
57 Mahunga Drive
Mangere Bridge, Auckland

Ph: 0800 800 823
sales@coopermedical.co.nz
www.coopermedical.co.nz

Diamond Scalpels

Handle Styles



President

Robust, all-titanium freehand diamond platform that pairs with paracentesis, phaco keratomes, crescents and pre-set series.

Spring-free mechanism, so there is no free play between hands and blade for ultimate precision.

Engineered to last with a lifetime handle warranty.



Elite II

Ultimate in controlled-depth precision for corneal relaxing incisions.

A premium Micrometer best paired with our Mastel Stealth Triamond diamond that has been geometrically engineered for better wound healing.



PHD II Micrometer

The "Swiss Army Diamond" of ophthalmic scalpels.

Simple, quick actuation by calibrated steps.

Includes full extension for hands-free surgery.

Diamond Scalpel Bestsellers

President Gimbel SuperStealth™ 2.0/2.5 Sharp Tip



SKU# 10.305.4033

- Gimbel SuperStealth™ provides surgeons with flexible incision profiles that are determined by how far diamond is advanced
- Sharp Enhancement Edges provide minimally invasive margins while maintaining ability to enlarge incision

Also Available:

10.305.4029 - 1.7/2.2mm Sharp

10.305.4031 - 1.7/2.2mm Bucci

10.305.4035 - 2.0/2.5mm Bucci

President UltraStealth 2.4mm Sharp Tip, Angled Mount



SKU# 10.305.4043

- Workhorse of any major cataract center
- Provides exceptional and consistent incision quality with industry leading UltraThin blade technology

Also Available:

10.305.4037 - 2.0mm Sharp

10.305.4039 - 2.9mm Sharp

10.305.4041 - 2.2mm Bucci

10.305.4045 - 2.4mm Bucci


President Meister Trifacet UltraThin Lance Thin, Straight Mount



SKU# 10.301.2013

- Parallel side edges allow for predictable incision size
- Utilizes Mastel's UltraThin Technology
- Truncated tip adds minimal resistance for paracentesis formation
- Available in straight or angled configurations

President Stealth Fine Triamond .3mm Tip, Angle Mount

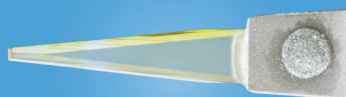

SKU# 10.302.2016

- Stealth's self-directing facet facilitates dive into anterior chamber during paracentesis formation
- Blade shape allows for variable incision width to match surgeon's preference
- Utilizes Mastel's UltraThin Technology to help facilitate wound closure

Also Available:

10.302.2018 President Stealth Enhanced Lance 60° Sharp Tip, Angled Mount

10.302.2019 President Stealth Enhanced Lance 60° Truncated Tip, Angled Mount

President Ahmed SuperCrescent 2.0 (M), 45° Angle Mount



SKU# 10.311.1087

- Crescent style 1.4mm blade with three facets offers one of the finest dissecting edges available
- Long cutting sides provide ability to enlarge incisions to surgeon's preference, and allows for better stability during flap creation

President Stealth Fine Triamond .3 Tip, Nichamin Classic 600 (.600 Preset)



SKU# 10.309.3016

- Popular blade configuration for LRI incisions due to the 7 deg blade tip angle which allows for more accurate incisional dimensions when compared to a standard spear shaped knife
- Polished footplate ensures a fixed and predictable depth for confidence inspiring outcomes

Also Available:

10.308.3016 President Stealth
Fine Triamond .3 Tip, Nichamin
Short .600 Preset

Elite II Fine Stealth Triamond .3 Tip, Corneal View



SKU# 12.322.1016

- Versatile diamond with long edges for maximum flexibility when making LRI's or groove incisions
- Handle style allows for variety of preset depths to ensure accuracy while consolidating what was historically several knives into one sleek package

Also Available:

12.321.1016 Insert Profile Used with Guide Ring

12.323.1016 Elite II Fine Stealth Triamond .3 Tip, Profile

12.323.2016 Elite II Fine Stealth Triamond .3 Tip, Profile, Compass

President EconoLance Standard, Straight Mount



SKU# 10.302.2001

- An economical sideport knife that allows for consistent incisional results
- Diamond shape allows for fixed incision width regardless of incision depth
- Available in both straight or angled configurations

Also Available:

10.302.2001 Angled Mount

President EconoLance Truncated (EconoMeister), Groover .550mm



SKU# 10.310.3086-55

- Exceptional low cost option for LRI or groove incisions
- Fixed depth allows for consistent incisions without need to adjust your knife

Also Available:

10.310.3086-25 Groover .250mm

10.310.3086-30 Groover .300mm

10.310.3086-35 Groover .350mm

10.310.3086-50 Groover .500mm



IMPORTANT: Mastel maintains strict quality standards for the dimensions of our diamond blades, and in certain circumstances a blade replacement will be required to ensure the consistent incisional quality of our knives.

Our Commitment



Mastel Surgical has proudly served the ophthalmology community for over 45 years. Our deep learnings from the operating room have enabled us to better understand surgeon needs and to develop superior products that support your dedication to delivering the best possible care to your patients.